

Capre MC Medical Chair

Patient comfort during demanding procedures like cytostatic and dialysis treatments



Product pictured with optional accessories.



Watch a video
<http://bit.ly/capretuoli>



Main features:

1. Modern solution to support the demanding treatments without the need of traditional hospital bed
2. An easy-to-use adjustments and pre-programmed positions
3. An eye-catching design and anatomically shaped upholstery enhance maximum patient comfort

The luxurious Capre multipurpose chair has been designed together with medical professionals to fulfil the demands of today's healthcare processes in demanding procedures such as cytostatic and dialysis treatments. The chair is used instead of the traditional hospital bed to improve patient comfort during treatments. The Capre chair features pleasantly contoured, modern upholstery, individual and easy to use adjustment possibilities (with up to five motors) and many optional built-in comforts, such as a tablet pc stand and reading light. All of this combines to improve patients' quality of life during treatments, whilst also improving the quality image of the healthcare facility.

The control system allows the adjustment of the chair height, back section, seat angle and leg section. Both the head section and the leg section also feature a length adjustment mechanism. The chair features a Trendelenburg position as well as three pre-programmed treatment positions (chair, beach chair and flat lying surface) and a user memory position in which the customer's own settings can be stored. The chair can be set, for example, to a lying down or resuscitation position by one press of a button. Standard multi-function armrests are height and angle adjustable as well as rotated. Optional sampling armrests are convenient for blood sampling, thanks to it wider, longer, and contoured shape of armrest upholstery. Motorised one-piece leg section with optional foot support gives extra assistance in dialysis treatments.

LOJER® For easy care


Technical data and features		Capre MC Multipurpose chair
Total length	186 - 213 cm	
Patient surface length	186 - 213 cm	
Total width	880 mm	
Patient surface width	620 mm	
Total weight	100 kg	
Height adjustment	60 - 85 cm	
Head section adjustment	10°, length adjustment 10 cm	
Back section adjustment	-12° - +90	
Back section length	600 mm	
Seat adjustment	14°	
Seat width	620 mm	
Leg section adjustment	-80° - 0°	
Arm support length	42 cm	
Arm support adjustment	Height adjustment 10 cm, angle adjustment -45° to +20° and rotation 0° - 180°	
Electric Trendelenburg	max. 14°	
Central locking castors	Ø75mm	
Other standard features	Motorised one-piece leg section USB charger Hand control with flexible arm or hand control with spiral cable Pre-programmed positions and one for user to decide White or grey cover for the base frame, under the seat and back section	
Upholstery	Anatomically contoured Future-upholstery, several colour options, for more info, please see Upholstery Product card	
Safe Working Load (SWL)	160 kg	
CE-marking	Yes	
Warranty	2 years	
Country of Manufacture	Finland	

Factory installed options:

- Electric foot control for height and back section adjustment
- Electric Trendelenburg; fixed foot control for Trendelenburg adjustments
- All foot controls (featuring both foot control for height and back section + fixed control for Trendelenburg)
- Battery for back-up
- Foot support for the standard motorised one-piece leg section

Other accessories:

- Sample arm support with anatomic upholstery (length 60 cm)
- Neck pillow
- Breathing hole
- Black tray
- Flexible tablet holder
- LED light
- IV-pole

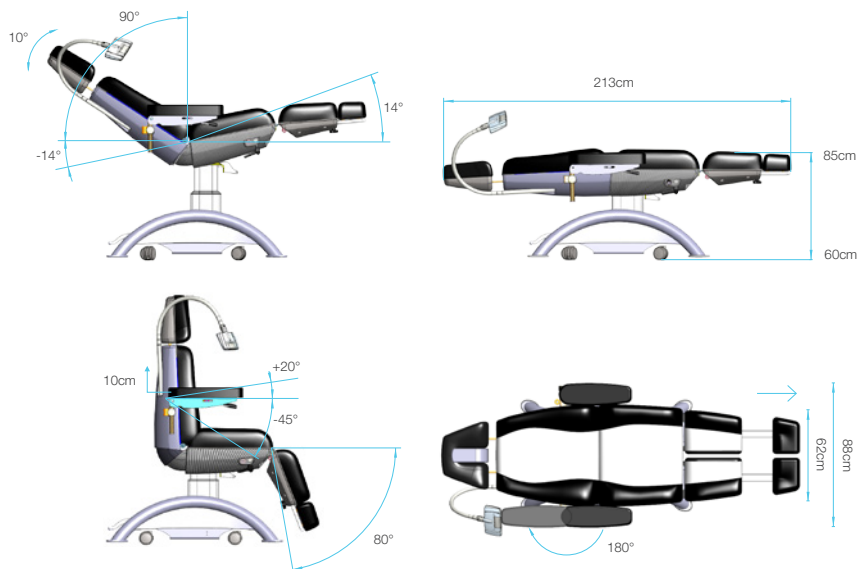

MADER
MEDICAL SYSTEMS

Kirchlergründe 6/19
A-6200 Jenbach

T: +43 5244 211 41 www.mader-medical.com
F: +43 5244 211 41-42 office@mader-medical.com