


Medical Recliner Chair 6801

Renewed and versatile hygiene patient chair



Product pictured with optional accessories.

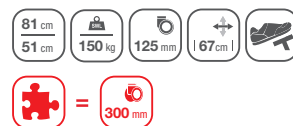

Main features:

1. Extensive range of hydraulic adjustments for the comfort of patients and care personnel
2. Easy to clean hygiene chair with seamless antimicrobial upholstery, hygiene handles and adjustment knobs
3. Comfortable to sit on - thanks to ergonomically shaped seat

The renewed hygiene patient chair 6801 is designed for medical examination, sampling and patient transfers. Special attention has been paid to the design, adjustability and ergonomics of the chair seat from both the nurse's and the patient's perspective. The chair has excellent mobility thanks to Ø125 mm castors. For more challenging conditions - such as outdoor use - you can choose optional Ø300mm front castors. Chair's upholstery is extremely durable and antimicrobial fabric that is also easy to clean, waterproof and has excellent disinfectant features.

The stepless back and leg section adjustments make the care work easier. The chair tilts into fully reclined position easily with the help of gas springs, and it can be lifted to the desired height hydraulically using the foot pump.

LOJER® For easy care

Technical data and features	
Total width	67 cm
Height adjustment (from the front edge of the seat)	51 - 81 cm
Seat width (between arm supports)	47 cm
Seat depth	46 cm
Back section adjustment	90°
Leg section adjustment	90°
Headrest height adjustment	9 cm
Castors	Ø125 mm locking castors including directional wheel
Safe working load (SWL)	150 kg
Weight	72 kg
Frame	Hygiene epoxy-coated steel frame
Upholstery	Seamless antimicrobial FUTURE upholstery; several color options
Other standard features	Folding foot rest, folding arm rests (5 positions with quick lock), ergonomically shaped seat, neck cushion, hygiene handles and adjusting knobs
CE marking	Yes
Warranty	2 years
Country of Manufacture	Finland

Factory installed options:

- Calf support length adjustment 130 mm

Other accessories:

- Arm rest for sampling
- Durable hygiene pillow for sampling 35x45cm (50045)
- IV-pole (60121AMI) and adapter
- Fixed Ø300 mm wheels in front



MADER
MEDICAL SYSTEMS

Kirchlergründe 6/19
A-6200 Jenbach

T: +43 5244 211 41 www.mader-medical.com
F: +43 5244 211 41-42 office@mader-medical.com